


# ENTERAL Catalogue



*Leaders in enteral delivery*

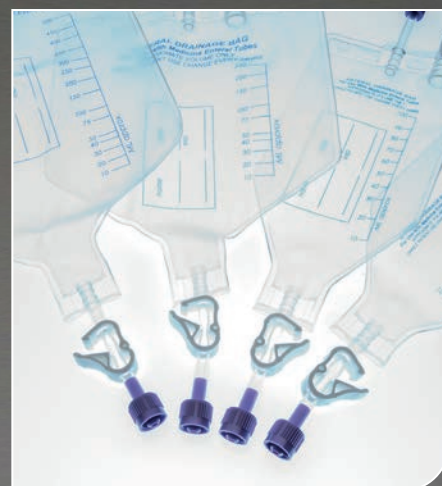

## Enteral Syringes

4-6

- ENFit Single Use Syringes
- ENFit Reusable Syringes
- Oral Tip Low Dose Syringes (ENFit compatible)
- Oral Tip Syringes
- Reusable Oral Tip Syringes

## Syringe Accessories

7-10

- ENFit Bottle Adapters
- ENFit Syringe Caps
- ENFit Straws
- ENFit Blunt Drawing Up Needles
- ENFit Spike Adapter
- ENFit Colostrum Collector
- ENFit Syringe Bolus Filter

## Nasogastric & Nasojejunal feeding tubes

11-15

- ENFit Polyurethane Nasogastric Feeding Tubes
- ENFit Polyurethane Nasogastric Feeding Tubes with Guidewire
- ENFit Polyurethane Y-Port Nasogastric Feeding Tubes with Guidewire
- ENFit Polyurethane Nasojejunal Feeding Tubes with Guidewire
- ENFit Feeding/ Aspiration Tubes

## Medicina Specialist Tubes

16-18

- ENFit Double Lumen Feeding / Aspiration Tubes
- ENFit Double Lumen Salem Tubes
- ENFit Double Lumen Replogle Tube
- Surfactant Kit

## ENFit Feeding Tube Accessories

19-26

- ENFit Drainage Bags
- ENFit Decompression System
- ENFit Gravity Sets
- ENFit Extension Sets
- Gastrostomy Bolus Extension Set
- ENFit 4 Way Stopcock
- pH Strips
- Nasogastric Nose Plaster
- Gastrostomy Tidy
- Ace Stopper
- ACElok Dressing
- Transitional Adaptors



Medicina has been supplying and innovating medical products for over 25 years.

We specialise in enteral feeding consumables including syringes, feeding tubes and accessories.



Medicina introduced the first safe enteral feeding system in 2005, the Reverse Luer system. In 2009 this system became the UK standard connector. The Reverse Luer system significantly improved patient safety and reduced the risk of accidental misconnections with other non-enteral patient access lines.

After seeing the systems advantages and safety improvements, the Reverse Luer system was adopted by additional countries until the global international design standard (ISO) connector was introduced.

Medicina successfully launched the ISO standard connector 'ENFit' (ISO80369-3) in 2015. The introduction of this standard marked a milestone in the global effort to improve patient safety. ENFit was designed to replace all other enteral connector systems used throughout the world to ensure one global solution. The ENFit standard establishes the requirements for enteral small-bore connectors making it difficult, if not impossible, for unrelated delivery systems to be connected. The ENFit system provides a "hospital-wide" solution for safer enteral access

because all components are compatible throughout the hospital. Medicina also offers a compatible and complimentary home feeding system ensuring a smooth transition from hospital to home.

Medicina was proud to become a part of the HMC group in 2017. HMC manufactures Medicina and own-brand medical devices in its state-of-the-art manufacturing facilities in Mirandola, Italy. Our shared vision is to provide practical, safe and reliable products to meet both clinician and patient needs. The combination of HMC's manufacturing expertise and Medicina's history of innovation within the enteral sphere, provides a great platform to provide leading enteral feeding devices, resulting in superior healthcare throughout the world.



*Leaders in enteral delivery*



# Enteral Syringes

The ENFit ISO Standard range of enteral syringes has been developed to deliver enteral feeds, medicines and flushes safely to neonatal, paediatric and adult patients.

This safety system provides a 'hospital wide' solution for safer enteral access because all the components are compatible with all other enteral feeding systems used within the hospital. A compatible home feeding system is also provided, ensuring a smooth transition from hospital to home. Oral tip syringes are also available to allow administration of feed/medications orally.



## ENFit Single Use Syringes

### Features

- Sterile
- Individually Packaged
- Purple
- ENFit compliant
- Clear measurement markings
- Low dose tip design to improve dose accuracy (1ml and 2.5ml)
- Three piece syringe (barrel, plunger and silicone gasket)
- DEHP/Latex free




Medicina code	Product Description	Box Qty
LPE01LD	1ml ENFit Low Dose Syringe	100
LPE25LD	2.5ml ENFit Low Dose Syringe	100
LPE05	5ml ENFit Syringe	100
LPE10	10ml ENFit Syringe	100
LPE20	20ml ENFit Syringe	80
LPE30	30ml ENFit Syringe	80
LPE60	60ml ENFit Syringe	60
LPE100	100ml ENFit Syringe	30





## ENFit Reusable Syringes

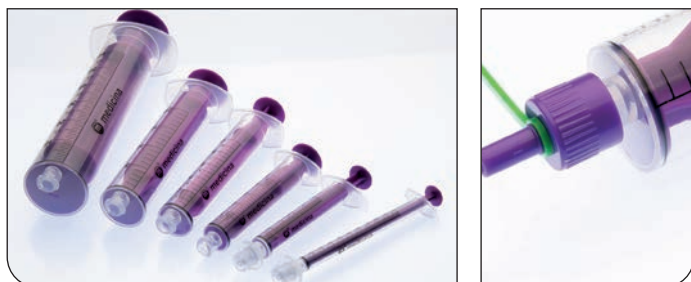

### Features

- Reusable/Washable
- Individually packaged
- Purple
- ENFit compliant
- Clear measurement markings
- Low dose tip design to improve dose accuracy (1ml and 2.5ml)
- Three piece syringe (barrel, plunger and silicone gasket)
- O-ring silicone gasket to prevent swelling when cleaning
- DEHP/Latex Free

Medicina code	Product Description	Box Qty
LHE01LD	1ml ENFit Reusable Low dose Syringe	100
LHE25LD	2.5ml ENFit Reusable Low Dose Syringe	100
LHE05	5ml ENFit Reusable Syringe	100
LHE10	10ml ENFit Reusable Syringe	100
LHE20	20ml ENFit Reusable Syringe	80
LHE60	60ml ENFit Reusable Syringe	60

## Oral Tip Low Dose Syringes (ENFit compatible)

The oral tip low dose syringe has been designed to give feed and medication orally, as well as being compatible with ENFit devices.



### Features

- Sterile
- Individually packaged
- Purple
- ENFit compatible
- Clear measurement markings
- Low dose tip design to improve dose accuracy
- Three piece syringe (barrel, plunger and silicone gasket)
- Smooth oral tip connector
- DEHP/Latex Free

Medicina code	Product Description	Box Qty
OT01LD	1ml Oral Tip Low-Dose Syringe (ENFit Compatible)	100
OT25LD	2.5ml Oral Tip Low Dose Syringe (ENFit Compatible)	100



## Oral Tip Syringes

### Features

- Sterile
- Individually packaged
- Purple
- Smooth oral tip connector
- Clear measurement markings
- Three piece syringe (barrel, plunger, and silicone gasket)
- DEHP/Latex Free



Medicina code	Product Description	Box Qty
OT005	0.5ml Oral Tip Syringe	100
OT01	1ml Oral Tip Syringe	100
OT25	2.5ml Oral Tip Syringe	100
OT05	5ml Oral Tip Syringe	100
OT10	10ml Oral Tip Syringe	100
OT20	20ml Oral Tip Syringe	80

## Reusable Oral Tip Syringes

### Features

- Reusable/Washable
- Individually packaged
- Purple
- Smooth oral tip connector
- Clear measurement markings
- Three piece syringe (Barrel, plunger and silicone gasket)
- DEHP/Latex Free



Medicina code	Product Description	Box Qty
OTH01	1ml Reusable Oral Tip Syringe	100
OTH25	2.5ml Reusable Oral Tip Syringe	100
OTH05	5ml Reusable Oral Tip Syringe	100
OTH10	10ml Reusable Oral Tip Syringe	100
OTH20	20ml Reusable Oral Tip Syringe	80



# Syringe Accessories

Medicina has put together a large range of accessories to help with the delivery of feeds and medicines using the ENFit syringe feeding system.



## ENFit Bottle Adapters

When drawing-up with an ENFit syringe, an ENFit Bottle Adapter allows for more accuracy as per ISO 80369-3 recommendation.

They are designed to fit tightly in the neck of the bottle allowing the tamperproof cap to be screwed on securely (excluding BA01L). ENFit Bottle Adapters will last the life of the bottle and they should be discarded together.



### Features

- Sterile
- Individually packaged (Excluding BAxxBL range)
- Purple
- ENFit compliant
- Available in various sizes to fit a range of bottles
- Stepped bottle adapter fits multiple bottles
- Minimises medicine waste
- Improves dose accuracy
- Tamper proof cap can still be utilised (excluding BA01L)
- DEHP/Latex Free

Medicina code	Product Description	Box Qty
BA01L	ENFit Stepped Bottle adapter	50
BA02L	ENFit Bottle Adapter 15mm-16.5mm	50
BA03L	ENFit Bottle Adapter 18mm – 19.5mm	50
BA04L	ENFit Bottle Adapter 25.5mm – 27mm	50
BA05L	ENFit Bottle Adapter 18.5mm -20.5mm	50
BA06L	ENFit Bottle Adapter 20mm – 21.5mm	50
BA07L	ENFit Bottle adapter 10.5mm -11mm	50
BA08L	ENFit Bottle Adapter 16.5mm – 18mm	50
BA01BL	ENFit Stepped Bottle Adapter (Pack)	50
BA02BL	ENFit Bottle Adapter 15mm – 16.5mm (Pack)	50
BA03BL	ENFit Bottle Adapter 18mm – 19.5mm (Pack)	50
BA04BL	ENFit Bottle Adapter 25.5mm 27mm (Pack)	40






## ENFit & Oral Tip Syringe Caps

ENFit/Oral Tip Syringe Caps allow feeds/medications to be prepared in advance.

ENFit/Oral Tip Syringes can be capped off until required. The syringe caps come packaged singularly or in packs of 8 to allow multiple feeds/medications to be prepared.

### Features

- Sterile
- Available singular or in packs of eight
- Purple
- ENFit compliant
- Oral Tip compliant (OTC range)
- Prevent syringes from leaking
- Connects with ENFit and ENFit low dose syringes
- Allows end users to prepare feeds/medications in advance
- DEHP/Latex Free

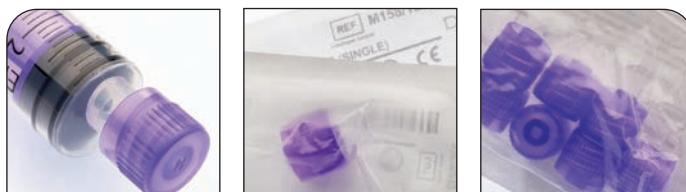

Medicina code	Product Description	Box Qty
M158/1LD	ENFit Syringe Cap (Single)	50
M158/8LD	ENFit Syringe Caps (50 packs of 8)	400
OTC01	Oral Tip Syringe Caps (Single)	50
OTC08	Oral Tip Syringe Caps (50 packs of 8)	400

## ENFit Straws

The Medicina ENFit Straws are available in various lengths and bore sizes, allowing you to draw up a wide range of different fluids.

### Features

- Sterile
- Individually packaged
- Purple connector
- ENFit compliant
- Oral Tip compliant (OTC range)
- Various lengths/bore sizes
- Minimises medicine waste
- Improves dose accuracy
- Filter straw is available for glass ampoules
- Larger straws prevent end users having to invert larger bottles
- Rectal straws for easy administration of rectal medications
- DEHP/Latex Free



Medicina code	Product Description	Box Qty
C005L	ENFit 5cm Medicine Straw	50
MF01L	ENFit Filter Straw	50
C010L	ENFit 10cm Straw	50
C014L	ENFit 14cm Straw	50
C020L	ENFit 20cm Straw	50
RC075L	ENFit 7.5cm Rectal Straw	50
RC150L	ENFit 15cm Rectal Straw	50





## ENFit Blunt Drawing Up Needles

When drawing up into an ENFit syringe from a vial with a rubber bung or a glass ampoule, you may need to use an ENFit Blunt fill needle or an ENFit Blunt fill filter needle.

The blunt needles are cut to a 45 degree angle, meaning they take 10 times the amount of force required to pierce human skin than a standard hypodermic needle. This allows you to access multi dose vials/ampoules whilst reducing the risk of needle stick injury.



### Features

- Sterile
- Individually packaged
- Purple connector
- ENFit compliant
- Reduces the risk of needle stick injury
- Complies with European Directive 2010/32/EU on sharps
- Reduces medicine waste
- Improves dose accuracy
- Sharp enough to access a rubber bung without coring
- 5 micron filter (MN02L)
- DEHP/Latex Free

Medicina code	Product Description	Box Qty
MN01L	ENFit Blunt Fill Needle	50
MN02L	ENFit Blunt Fill Filter Needle	50

## ENFit Spike Adapter

The ENFit Spike Adapter allows multiple ENFit syringe access to feed bags.



### Features

- Sterile
- Individually packaged
- Purple connectors
- ENFit compliant
- ENPlus connection
- DEHP/Latex Free

Medicina code	Product Description	Box Qty
SA01L	ENFit Spike Adapter	50



## ENFit Colostrum Collector

When the ENFit standard was introduced, especially with the low dose amendment, it became almost impossible to express colostrum into smaller syringes due to the spike in the low dose connector.

When you use the colostrum collector in conjunction with OTLD syringes (page 5) you have the ideal solution.

### Features

- Sterile
- Individually packaged
- Purple connector
- ENFit compliant
- Easy to collect colostrum
- Comfortable for patient
- DEHP/Latex Free



Medicina code

Product Description

Box Qty

TA10

ENFit Colostrum Collection Device

50

## ENFit Syringe Bolus Filter Cap

The Bolus filter cap protects feed from contamination. It can be used in conjunction with 20/30/60ml ENFit syringes.

### Features

- Sterile
- Individually packaged
- Purple
- ENFit compliant
- Protects feed from contamination
- DEHP/Latex free



Medicina code

Product Description

Box Qty

BA10L

ENFit Syringe Bolus Filter

50



# Nasogastric & Nasojejunal feeding tubes

Medicina has a comprehensive range of nasogastric and nasojejunal tubes to feed a wide variety of patients from neonates to adults both in hospital and at home.

## Fully NPSA compliant

In March 2011 the NPSA UK issued a Patient Safety Alert which updated the NPSA Alert from March 2005. These alerts are very clear in the physical characteristics needed for a tube that is going to be used for feeding, flushing or for delivering medication. All Medicina nasogastric and nasojejunal tubes are fully compliant with both alerts.


## ENFit Polyurethane Nasogastric Feeding Tubes



Medicina code	Product Description	Box Qty
LG4/40	4FR x 40cm	50
LG4/50	4FR x 50cm	50
LG4/80	4FR x 80cm	50
LG5/40	5FR x 40cm	50
LG5/50	5FR x 50cm	50
LG5/80	5FR x 80cm	50
LG6/40	6FR x 40cm	50
LG6/50	6FR x 50cm	50
LG6/80	6FR x 80cm	50
LG8/40	8FR x 40cm	50
LG8/50	8FR x 50cm	50
LG8/80	8FR x 80cm	50
LG8/120	8FR x 120cm	50
LG10/80	10FR x 80cm	50
LG10/120	10FR x 120cm	50
LG12/80	12FR x 80cm	50
LG12/120	12FR x 120cm	50
LG14/80	14FR x 80cm	50
LG6/50W	6FR x 50cm (weighted)	50
LG6/80W	6FR x 80cm (weighted)	50
LG8/50W	8FR x 50cm (weighted)	50
LG8/80W	8FR x 80cm (weighted)	50
LG10/80W	10FR x 80cm (weighted)	50
LG12/80W	12FR x 80cm (weighted)	50

## Features

- Sterile
- ENFit compliant
- Polyurethane
- Radio-opaque stripe
- Centimetre markings throughout the length of the tube
- Two side eye holes
- Colour coded strap to signify french size of tube
- Tethered cap
- Rounded tip to aid insertion
- Weighted tubes available
- Up to 90 days use

 Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex Free

UpTo  
**90**

Additional size tubes are available upon request.



# ENFit Polyurethane Nasogastric Feeding Tubes with Guidewire



## Features

- Sterile
- ENFit compliant
- Polyurethane
- Fully radiopaque
- Barium sulphate tip – highlights end of tube clearly under X-Ray
- Centimetre markings throughout the length of the tube
- Two side eye holes
- Colour coded strap to identify French size of tube
- Tethered cap
- Rounded tip to aid insertion
- Guidewire to aid insertion (ENFit connection)
- Includes documentation label
- Weighted tubes available
- Up to 90 days use



Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex Free



Medicina code	Product Description	Box Qty
<b>Nasogastric Tubes (with guidewire)</b>		
NGP6/55L	6FR x 55cm	25
NGP6/75L	6FR x 75cm	25
NGP6/85L	6FR x 85cm	25
NGP8/55L	8FR x 55cm	25
NGP8/75L	8FR x 75cm	25
NGP8/85L	8FR x 85cm	25
NGP8/120L	8FR x 120cm	25
NGP10/85L	10FR x 85cm	25
NGP10/120L	10FR x 120cm	25
NGP12/100L	12FR x 100cm	25
NGP12/120L	12FR x 120cm	25
NGP14/120L	14FR x 120cm	25
<b>Nasogastric Tubes Weighted (with guidewire)</b>		
NGP8/85LW	8FR x 85cm	25
NGP10/85LW	10FR x 85cm	25

UpTo  
**90**





# ENFit Polyurethane Y-Port Nasogastric Feeding Tubes with Guidewire

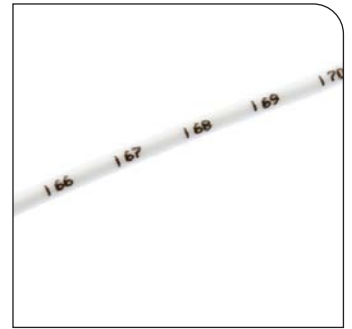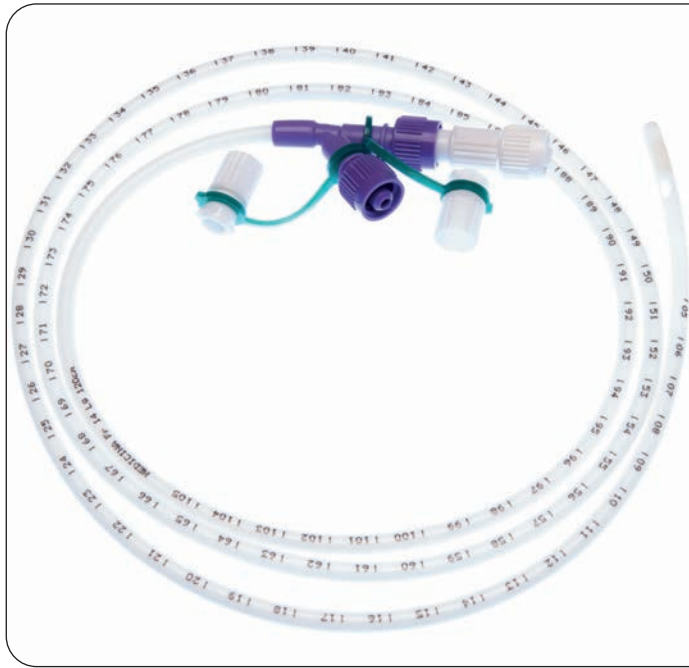

Medicina code	Product Description	Box Qty
<b>Y-Port Nasogastric Tubes (with guidewire)</b>		
NGP5/55LY	5FR x 55cm	25
NGP5/75LY	5FR x 75cm	25
NGP6/55LY	6FR x 55cm	25
NGP6/75LY	6FR x 75cm	25
NGP8/55LY	8FR x 55cm	25
NGP8/75LY	8FR x 75cm	25
NGP8/85LY	8FR x 85cm	25
NGP8/10LY	8FR x 100cm	25
NGP8/11LY	8FR x 110cm	25
NGP8/12LY	8FR x 120cm	25
NGP10/85LY	10FR x 85cm	25
NGP10/10LY	10FR x 100cm	25
NGP10/11LY	10FR x 110cm	25
NGP10/12LY	10FR x 120cm	25
NGP12/85LY	12FR x 85cm	25
NGP12/10LY	12FR x 100cm	25
NGP12/11LY	12FR x 110cm	25
NGP12/12LY	12FR x 120cm	25
NGP14/85LY	14FR x 85cm	25
NGP14/10LY	14FR x 100cm	25
NGP14/12LY	14FR x 120cm	25

## Features

- Sterile
- Y-Port ENFit compliant
- Polyurethane
- Fully radiopaque
- Barium sulphate tip – highlights end of tube clearly under X-Ray
- Centimetre markings throughout the length of the tube
- Two side eye holes
- Colour coded strap to identify French size of tube
- Tethered cap
- Rounded tip to aid insertion
- Guidewire to aid insertion (ENFit connection)
- Includes documentation label
- Up to 90 days use



Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex Free

UpTo  
**90**



# ENFit Polyurethane Nasojejunal Feeding Tubes with Guidewire

The nasojejunal tubes have an open end feeding port that allows them to be passed over an endoscopically or radiologically placed guidewire.



## Features

- Sterile
- ENFit compliant
- Polyurethane
- Fully radiopaque
- Centimetre markings throughout the length of the tube
- One side eye hole
- Colour coded strap to identify French size of tube
- Tethered cap
- Open tip (Excluding weighted tubes)
- Guidewire to aid insertion (ENFit connector)
- Includes documentation label
- Weighted tubes available
- Up to 90 days use



Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex Free

Medicina code	Product Description	Box Qty
NJP6/120L	6FR x 120cm	10
NJP8/120L	8FR x 120cm	10
NJP8/130L	8FR x 130cm	10
NJP10/130L	10FR x 130cm	10
NJP10/130LW	10FR x 130cm (Weighted)	10
NJP12/130L	12FR x 130cm	10
NJP12/130LW	12FR x 130cm (Weighted)	10
NJP14/130L	14FR x 130cm	10
NJP14/130LW	14FR x 130cm (Weighted)	10

Up To



90



# ENFit Feeding/Aspiration Tubes


These tubes are primarily designed to aspirate stomach contents in adults and children.

They can also be used for continuous or intermittent nasogastric enteral feeding, administering medication and flushing.

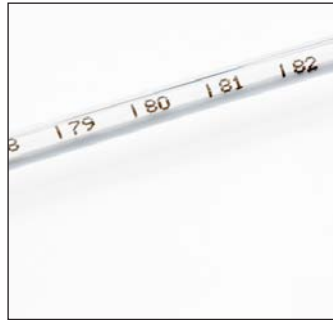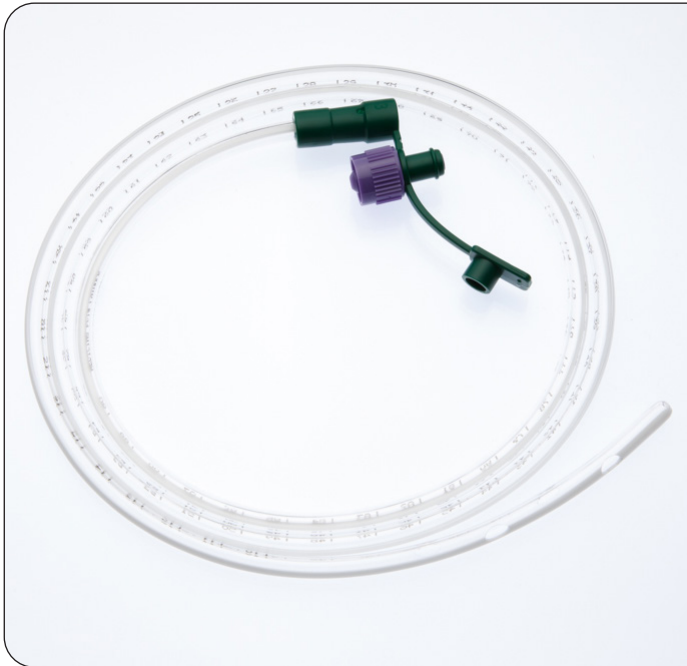

## Features

- Sterile
- ENFit compliant
- Polyurethane
- ENFit and funnel connector
- Radio-opaque stripe
- Centimetre markings throughout the length of the tube
- Four large side eye holes
- Colour coded strap to signify french size of tube
- Tethered cap
- Rounded tip to aid insertion
- Includes documentation label
- Weighted tubes available
- Tube including guidewire available
- Up to 90 days use

 Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex Free

Medicina code	Product Description	Box Qty
RT6/80L	6FR x 80cm	25
RT8/80L	8FR x 80cm	25
RT10/100L	10FR x 100cm	25
RT12/100L	12FR x 100cm	25
RT14/100L	14FR x 100cm	25
RT16/100L	16FR x 100cm	25
RT18/100L	18FR x 100cm	25
RT20/100L	20FR x 100cm	25

Up To  
**90**



# Medicina Specialist Tubes




## ENFit Double Lumen Feeding/Aspiration Tubes

The ENFit double lumen nasojejunal feeding and gastric aspiration tube is used to feed into the small bowel and simultaneously aspirate the stomach. There are two versions of the tube depending on the method of placement.

The EF2L has an internal guidewire for surgical or endoscopic placement.

The EF2EL has an external guidewire and an open tip for an "Over the Wire" placement technique either radiologically or by using the guidewire down the centre of the endoscope.

The aspiration lumen of these tubes may be used to feed the patient as they recover. Care must be taken changing the labels on the tubes.

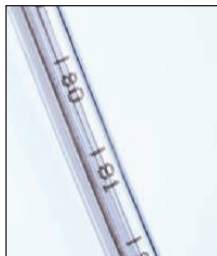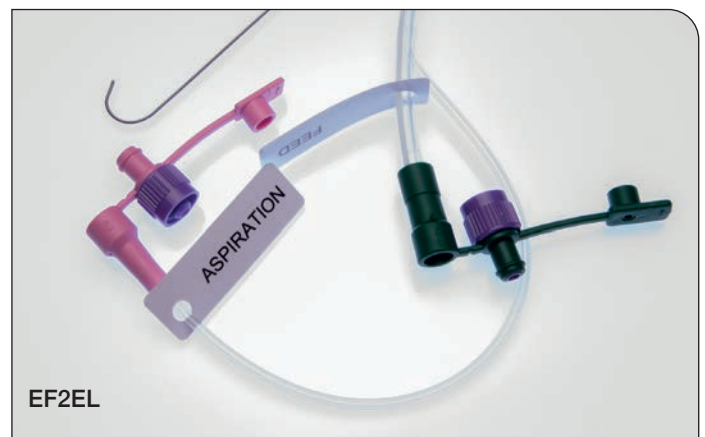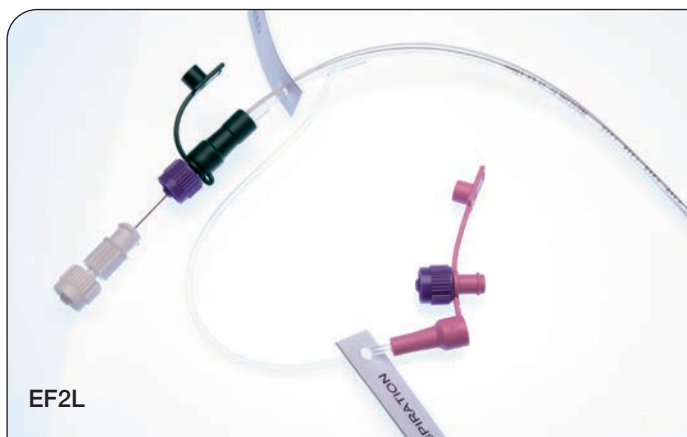

### Features

- Sterile
- ENFit compliant
- Double lumen tube
- Polyurethane
- Radio-opaque stripe
- Centimetre markings throughout the length of the tube
- Internal or external guidewire depending on method of placement
- Rounded tip or open ended tip depending on method of placement
- Up to 90 days use

 Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex Free

Medicina code	Product Description	Box Qty
EF2L	ENFit Double Lumen Feeding / Aspiration Tube with internal guidewire and closed tip	5
EF2EL	ENFit Double Lumen Feeding / Aspiration Tube with external guidewire and an open tip	5


UpTo  
**90**



# ENFit Double Lumen Salem Tubes


The “SALEM” type tubes are designed for simultaneous gastric aspiration and ventilation.

The double lumen limits depression of the stomach and protects the stomach lining. They can also be used for continuous or intermittent nasogastric enteral feeding, administering medication and flushing.

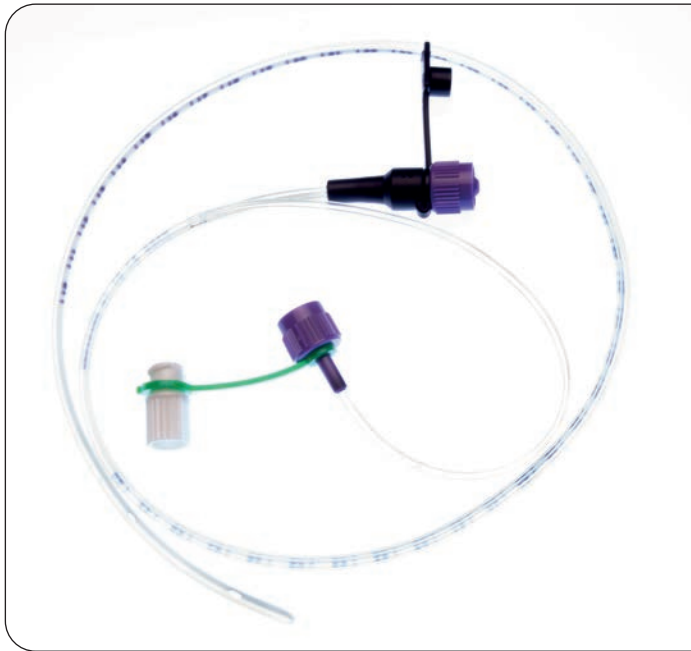

## Features

- Sterile
- ENFit Compliant
- Double lumen for Aspiration/ Ventilation
- Polyurethane
- Full length radiopaque stripe
- Centimetre markings throughout the length of the tube
- Four large side eye holes
- Rounded tip to aid insertion
- Colour coded to signify french size of tube
- Up to 90 days use

 Hospital guidelines must be followed for nasal erosion and ventilated patients.

- DEHP/Latex free

Medicina code	Product Description	Box Qty
25010ML	10FR x 60cm ENFit Salem Tube	10
25012ML	12FR x 120cm ENFit Salem Tube	10
25014ML	14FR x 120cm ENFit Salem Tube	10
25016ML	16FR x 120cm ENFit Salem Tube	10

UpTo  
**90**



# ENFit Double Lumen Replogle Tube

ENFit Double Lumen Replogle Tubes are used in Neonatal units to aspirate and irrigate fluids in conditions such as oesophageal atresia.

## Features

- Sterile
- Funnel connector for aspiration (black)
- ENFit compatible irrigation port (white)
- Polyurethane
- Double lumen allowing simultaneous irrigation and aspiration
- Radio-opaque stripe
- Tethered cap
- Centimetre markings throughout the length of the tube
- Rounded tip to aid insertion
- DEHP/Latex Free



Not for enteral administration



Medicina code	Product Description	Box Qty
NF02L	ENFit 10FR Double Lumen Replogle Tube	10

# Surfactant Kit

Sterile kit to easily and safely draw up and administer surfactant. The surfactant kit contains all the equipment necessary for this procedure.

## Features

- Sterile
- Dedicated surfactant kit connection
- Coloured blue to differentiate from feeding devices
- Complete kit.
  - 4FR 20cm Tube
  - 5ml syringe
  - Blunt Needle or Needle Free device
- Complies with EU directive on Sharps Safety
- DEHP/Latex Free



Not for enteral administration



Medicina code	Product Description	Box Qty
SF01	Surfactant Kit	25
SF02	Needle Free Surfactant Kit	25



# ENFit Feeding Tube Accessories

## ENFit Drainage Bags

Enteral drainage bags connect directly to feeding tubes or to ENFit extension sets (connected to feeding devices) to allow collection of drained fluids.



Medicina code	Product Description	Box Qty
MED50L	50ml ENFit Drainage Bag	25
MED100L	100ml ENFit Drainage Bag	25
MED250L	250ml ENFit Drainage Bag	25
MED500L	500ml ENFit Drainage Bag	25
MED1000L	1000ml ENFit Drainage Bag	25
MED51L	50ml ENFit Drainage Bag Including 75cm extension tubing	25
MED101L	100ml ENFit Drainage Bag Including 75cm extension tubing	25
MED251L	250ml ENFit Drainage Bag Including 75cm extension tubing	25
MED501L	500ml ENFit Drainage Bag Including 150cm extension tubing	25
MED1001L	1000ml ENFit Drainage Bag Including 150cm extension tubing	25

### Features


- Sterile
- Individually packaged
- Purple connector
- ENFit connector
- Extract samples of gastric fluid for analysis
- Clamp to prevent leaking
- Antibacterial filter allows gas to escape so the bag collects 100% stomach contents
- Five sizes are available
- DEHP/Latex free



# ENFit Decompression System

The ENFit Gastric Decompression System removes gas from a patient's stomach, eliminating the pain and discomfort associated with fluid and gas retention.

Patients can digest enteral formula at their own rate and reach their calorific goals.

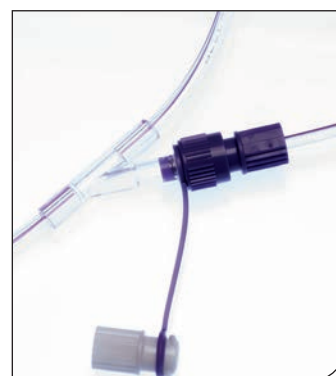

## Features

- Sterile
- Individually packaged
- Purple connector
- ENFit connector
- Relieves gastric distention (bloating of the stomach)
- Eliminates the pain and discomfort associated with fluid and gas retention
- Allows patients to digest enteral formula at their own rate and reach their enteral calorific goals
- Available in two sizes (Children & Adults)
- DEHP/Latex free

Medicina code	Product Description	Box Qty
EV01L	250ml ENFit Gastric Decompression System	30
EV02L	500ml ENFit Gastric Decompression System	30





# ENFit Gravity Sets

Medicina have a range of ENFit Gravity Feed Sets including a Universal Bottle Adapter or ENFit Cross Spike.



## Features


- Sterile
- Individually packaged
- Purple connectors
- ENFit connector
- Available with Universal Bottle Adapter or ENFit Cross Spike
- Connects to ENFit feeding tubes
- Drip Chamber
- ENFit Y Site
- Roller Clamp
- The MGS12L and MGS13L include a detachable Funnel Adapter. The detachable Funnel Adapter can become a part of a gastrostomy feeding device making it ENFit compatible
- DEHP/Latex free

Medicina code	Product Description	Box Qty
MGS10L	ENFit Gravity Set with Universal Bottle Adapter	30
MGS11L	ENFit Gravity Set with ENFit Cross Spike	30
MGS12L	ENFit Gravity Set with Universal Bottle Adapter & Detachable Funnel Adapter	30
MGS13L	ENFit Gravity Set with ENFit Cross Spike & Detachable Funnel Adapter	30



# ENFit Extension Sets

Several long and short ENFit extension sets are available (with and without Y-Sites) to connect an ENFit device/pump to a feeding tube.

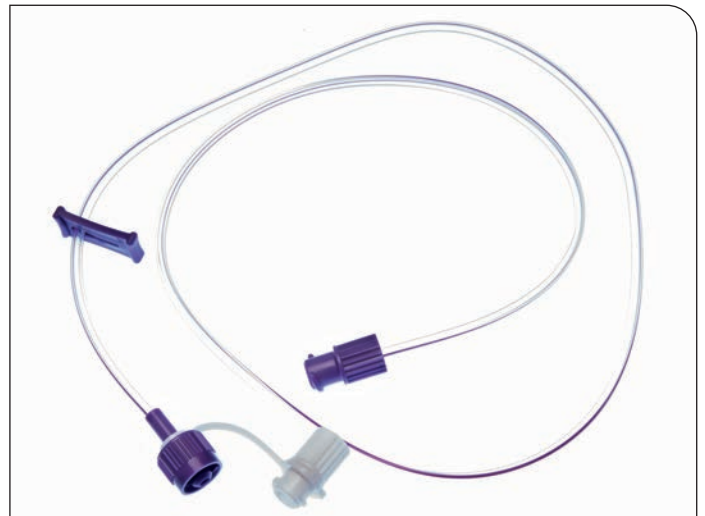

## Features


- **Sterile**
- **Individually packaged**
- **Purple connectors**
- **ENFit connectors**
- **Extension sets provide longer length tubing to connect to a feeding tube**
- **Wide bore extension sets are available for high density feeds**
- **Y-site option for administration of medication**
- **DEHP/ Latex Free**

Medicina code	Product Description	Box Qty
EP015L	15cm Long Extension Set ENFit	50
EP015LYW	15cm Long Y Site Extension Set ENFit Wide Bore	50
EP075L	75cm Long Extension Set ENFit	50
EP075LY	75cm Long Y Site Extension Set ENFit	50
EP075LW	75cm Long Extension Set ENFit Wide Bore	50
EP075LYW	75cm Long Y Site Extension Set ENFit Wide Bore	50
EP150L	150cm Long Extension Set ENFit	50
EP150LY	150cm Long Y Site Extension Set ENFit	50



## Gastrostomy Bolus Extension Set

The Gastrostomy Bolus Extension Set aids the process of feeding down a gastrostomy tube. It is used when the rate of delivery needs to be controlled.

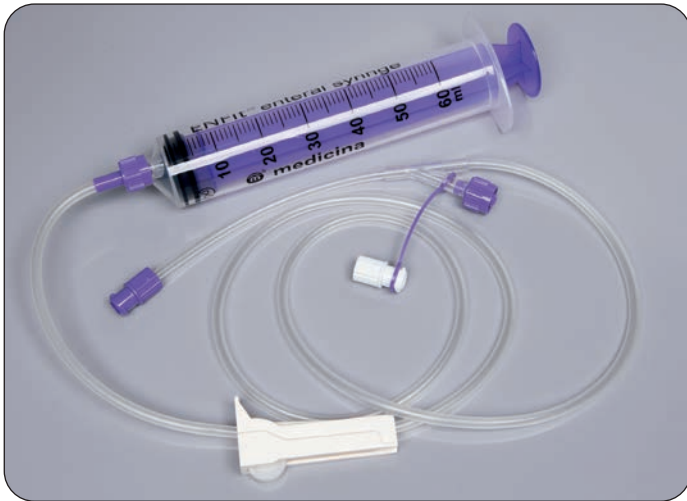

### Features

- Sterile
- Individually packaged
- Purple connectors
- ENFit connectors
- Roller clamp to control speed of delivery
- Y-site for administration of medication
- 60ml ENFit syringe included
- DEHP/ Latex free

Medicina code	Product Description	Box Qty
MGS02L	Bolus or Gravity Gastrostomy Extension Set	30

## ENFit 4 Way Stopcock

The ENFit 4-way stopcock has a female connector and two male connectors.

The device is used to provide an alternative route to administer enteral medication down a feeding tube. It can also be used to mix enteral medication and saline solution, in line with pharmacy advice and manufacturer's instructions.



### Features

- Sterile
- Individually packaged
- Purple connectors
- ENFit connectors
- Provides an alternative route to administer medication down a feeding tube
- Can be used to mix enteral medication and saline solution
- DEHP/ Latex free

Medicina code	Product Description	Box Qty
SA03L	ENFit 4 - way stopcock (one female port & two male)	50



## pH Strips

Under the NPSA guidelines, pH testing is the first line test for tube position.

### Features

- CE marked for testing human gastric aspirate
- Buffered against water
- The Medicina pH Indicator strip acts as its own colour chart, you do not need access to a box/ chart
- Available in boxes of 100 or 200



Medicina code	Product Description	Box Qty
IN01	pH Indicator Strips range 1-12	200
IN02	pH Indicator Strips range 1-12	100

## Nasogastric Nose Plaster

The Nasogastric Nose Plaster is a fixation device for Nasogastric and Nasojejunal tubes.

### Features

- Secures NG/ NJ tubes at nose or below the cheek bone
- Comfortable for patient
- Available in two sizes (Children & Adults)



Medicina code	Product Description	Box Qty
NP01	Adult Nose Plaster	100
NP02	Paediatric Nose Plaster	100

## Gastrostomy Tidy

The gastrostomy tidy is a device designed for patients who are body conscious or who pull on their gastrostomy tubes. The device allows a gastrostomy tube to be kept tidy and secure in a specially made pouch.

### Features

- Soft and comfortable to wear
- Discreet when worn underneath clothing
- Adjustable waist belt
- Machine washable
- Available in child and adult sizes



Medicina code	Product Description	Box Qty
GT01	Medium Gastrostomy Tidy (75cm - 114cm)	1
GT02	Large Gastrostomy Tidy (117cm - 152.5cm)	1



## Ace Stopper

Medicina ACE stoppers are small silicone plugs that are used to seal stoma sites at various parts of the body.

Ace stoppers can be used when a gastrostomy device fails or needs replacing. They will seal and maintain the patency of the stoma. A lot of healthcare professionals choose Ace Stoppers as a back-up for all gastrostomy fed patients.




### Features

- Sterile
- Individually Packaged
- Soft and comfortable
- Various lengths and shaft diameters depending on the size of the stoma
- Prevents leakage
- Seals the Stoma
- One ACElok dressing included with each ACE Stopper
- Available on prescription (UK)
- DEHP /Latex Free

Medicina code	Product Description	Box Qty
AP8/15	8FR - 15mm long Ace Stopper	1
AP8/30	8FR - 30mm long Ace Stopper	1
AP8/60	8FR - 60mm long Ace Stopper	1
AP10/15	10FR - 15mm long Ace Stopper	1
AP10/30	10FR - 30mm long Ace Stopper	1
AP10/60	10FR - 60mm long Ace Stopper	1
AP12/30	12FR - 30mm long Ace Stopper	1
AP12/60	12FR - 60mm long Ace Stopper	1
AP12/100	12FR - 100mm long Ace Stopper	1
AP14/30	14FR - 30mm long Ace Stopper	1
AP14/60	14FR - 60mm long Ace Stopper	1
AP14/100	14FR - 100mm long Ace Stopper	1

## ACElok Dressing

The ACElok dressing has been specifically designed for placement over the ACE stopper to provide external fixation and cushioning.

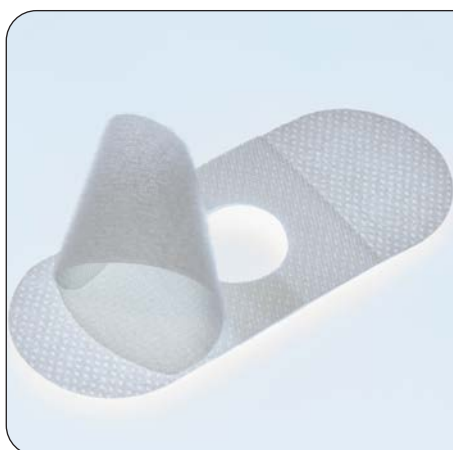

### Features

- Individually Packaged
- Used in combination with ACE Stoppers
- Allows easy access to a stoma site without having to remove dressing
- Lasts 3-4 days
- Available on prescription (UK)

Medicina code	Product Description	Box Qty
SD03	ACElok Dressing	10





## Transitional Adapters

The change to the globally recognised ENFit connector will mean that users will have to run down stock of “old” products before moving to the ENFit system.

All Medicina adapters are designed to be connected to the end of feeding tubes and can be used for up to 30 days.



Medicina code	Product Description	Box Qty
TA01	Transitional Adapter for Reverse Luer to ENFit Syringe	50
TA02	Transitional Adapter for ENFit Tube to Reverse Luer Syringe	50
TA03	Transitional Adapter for Large Cone Tube to ENFit Syringe/Set	50
TA04	Transitional Adapter for ENFit Tube to Large Cone Syringe	50
TA05	Transitional Adapter for ENFit Tube to Enlock Giving Set/Syringe	50
TA06	Transitional Adapter for ENFit Tube to Oral Tip Syringe/Set	50
ML01L	Transitional Adapter for Oral Tip Tube to ENFit Syringe/Set	50
LK01	Transitional Adapter for Funnel Tube to ENFit Syringe	50



Visit our website at [www.medicina.co.uk](http://www.medicina.co.uk)

Brochures and resources are available to download.

medicina.co.uk

SEARCH +44 (0)1204 695050

**m** medicina  
part of HMC Group  
Leaders in enteral delivery

Corporate Recruitment News Contact Us

ENFIT ORAL TIP SAFETY & IV STOMA CARE SILO BAGS GASTROSTOMY SURFACTANT KIT RESOURCES SALES CONTACTS

# ENFit<sup>®</sup>

The new International Enteral Safety Standard

**m** medicina  
**ENFit**

Medicina syringe feeding system

The worldwide introduction of the new ENFit connector will reduce the risk of enteral tube feeding misconnections and improve patient safety throughout hospitals around the globe.

**m** medicina  
**ENFit**

Medicina nasogastric feeding tubes

Medicina nasogastric feeding tubes are fully NPFA compliant and have all of the physical characteristics required for feeding, flushing and delivering medication.

The worldwide introduction of the new ENFit connector will reduce the risk of enteral feeding tube misconnections and improve patient safety throughout hospitals and around the globe.

**m** medicina  
**IVsyringes**

Medicina IV syringes

Cost saving solutions without compromising safety and efficacy.

**m** medicina  
**STOMACare**

Medicina developed this range of stoma care products specifically to meet new demands in the care of patients with gastrointestinal and urological stomas.

**m** medicina  
**Silo**

Medicina Silo  
An innovative surgical solution for infants with Gastrochisis.





Medicina has been supplying  
and innovating medical products  
for over 25 years.



*Leaders in enteral delivery*

Units 2-4, Rivington View Business Park,  
Station Road, Blackrod, Bolton BL6 5BN, UK  
Telephone: +44(0)1204 695050

[www.medicina.co.uk](http://www.medicina.co.uk)



[www.facebook.com/MedicinaUK](https://www.facebook.com/MedicinaUK)



[@medicina\\_ltd](https://twitter.com/medicina_ltd)

