

DIMAR

PRODUCT LINE

Intensive Care

DimAir Line



ECO OBLÒ HELMETS
with integrated safety valve

ECO OBLÒ HELMETS


with integrated safety valve

Single-patient helmet for non-invasive respiratory support in CPAP and BI-PAP with high compliance, complete with integrated two-way safety valve.

- Quick and direct access to the patient via removable screw-type airlock.
- The oblò facilitates all nursing, feeding and weaning procedures without having to remove the device.
- The oblò has a TWO-WAY anti-suffocation valve that allows the patient to exchange air with the outside (inlet and outlet) in case of failure to supply ventilation gas. The device is automatically activated when the pressure inside the helmet is less than 2.5 cmH₂O approximately.
- Available the version with valve equipped with a filter membrane with efficiency >99% , suitable to limit the spread of micro-organisms in the environment.

TECHNICAL FEATURES:

- Single patient
- Air-tight oblò with two-way antisuffocation valve with or without filter material
- Made of biocompatible plastic material
- Soft and flexible armpit straps with protection in hygienic antibacterial material CoolMax®
- Connectors:
 - 4 connectors ISO on the rigid ring
 - 2 lateral connectors, 22M, 22F ISO
- 2 retaining loops on the top of the helmet for an alternative/additional fastening
- Latex and PHT Free

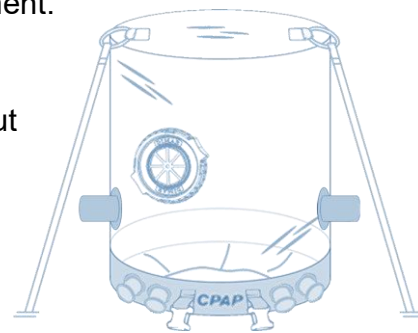

The use of the alternative/additional "eyelet" attachment prevents decubitus to the armpits



ADVANTAGES

- Effective, practical and comfortable
- Rapid patient access
- Possibility to perform nursing, feeding and weaning without removing the helmet
- TWO-WAY ANTI-SUFFOCATION valve for patient safety
- The version with filter material integrated in the valve is suitable for limiting the spread of micro-organisms in the environment

REFERENCES:

- In single pack
- Multiple pack of 5
- Possibility to create models, versions and custom packaging on specific request.

*+/- dimensional tolerances excluding collar and patient head

| CODE two-way valve WITH filter | CODE two-way valve WITHOUT filter | SIZE | Neck circumference (cm) | RDM | CND/EMDN |
|--------------------------------------|---|------|-------------------------------|---------|-------------|
| 500/6168 | 500/6061 | XXS | 23 ÷ 27 | 2683503 | R0301010401 |
| 500/6169 | 500/6062 | XS | 27 ÷ 31 | | |
| 500/6170 | 500/6063 | S | 31 ÷ 35 | | |
| 500/6171 | 500/6064 | M | 35 ÷ 39 | | |
| 500/6172 | 500/6065 | L | 39 ÷ 43 | | |
| 500/6173 | 500/6066 | XL | 43 ÷ 47 | | |
| 500/6174 | 500/6067 | XXL | 47 ÷ 51 | | |
| 500/6175 | 500/6068 | 3XL | ≥ 52 | | |



HMC Premedical S.p.A

Via Bosco 1/3 – 41037 Mirandola – MO – I
Tel: +39 0535 22704 Fax: +39 0535 609546

