

DIMAR

PRODUCT LINE

Ambulance
Emergency
& Hospital Departments

DimAir Line




***SYSTEM FOR
POSITIVE PRESSURE***

Positive pressure system formed by single tube Venturi generator, “T” piece assembled with manometer and mouthpiece. Single patient Venturi system, able to generate a FiO2 of approx. 0,35, equipped with 180 cm long O2 hose and connector.

TECHNICAL FEATURES:

- 1. Mouthpiece
- 2. T Piece
- 3. Single tube venturi generator
- 4. Manometer
- 5. Oxygen tube 180 cm with connector

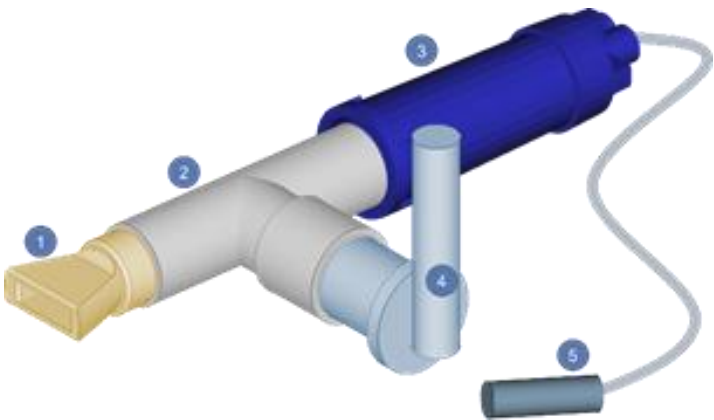

ADVANTAGES:

- Device intended for the treatment of respiratory insufficiencies in a non invasive way
- Expand lungs by increasing the capacity of functional residues (FRC)
- Treatment of atelectasia



- In single pack
- Multiple pack of 10

The system can be used with the EASY VENT VENTURI AEROSOL KIT

DESCRIPTION	SPU	CND/EMDN	RDM
SYSTEM FOR POSITIVE PRESSURE	700/3424	R03010180	2683488

DESCRIPTION	SPU	CND/EMDN	RDM
AEROSOL KIT	700/5066	N.A.	142843

N.A.
For assembled Art. 22 - MDR 2017/745/UE the CND is not provided, which is indicated for the individual products components of the KIT and contained within the configuration, as shown in the individual data sheets.